

# wellcare



## AIR COMPRESSION FOOT MASSAGER

### **USER MANUAL**

**Model Number FE-7203**



# Contents

03	<b>Important Safety Instructions</b>
07	<b>Product Overview</b>
08	<b>Operation</b>
08	<b>Troubleshooting</b>
09	<b>Specifications</b>

# Important Safety Instructions

## **READ CAREFULLY AND KEEP FOR FUTURE REFERENCE**

Read this manual thoroughly before first use, even if you are familiar with this type of product. The safety precautions enclosed herein reduce the risk of fire, electric shock and injury when correctly adhered to. Make sure you understand all instructions and warnings. Keep the manual in a safe place for future reference, along with any warranty information, your purchase receipt and the original packaging carton. If you sell or transfer ownership of this product, pass on these instructions to the new owner of the product.

**Always follow basic safety precautions and accident prevention measures, including the following:**

### **Usage conditions and restrictions**

- This product is intended for private use by adults only. It is not suitable for commercial or trade use. Do not use it for anything other than its intended purpose, and only use it as described in this manual.
- **ALWAYS** unplug this appliance from the electrical outlet immediately after using and before cleaning.
- **DO NOT** reach for product that has fallen into water. Unplug immediately.
- **DO NOT** place or store appliance where it can fall or be pulled into a tub or sink.
- **DO NOT** place in or drop into water or other liquid.

## Important Safety Instructions (Cont.)

### **WARNING: TO REDUCE THE RISK OF BURNS, FIRE, ELECTRIC SHOCK OR INJURY TO PERSONS:**

- An appliance should never be left unattended when plugged in. Unplug from outlet when not in use and before putting on or taking off parts.
- Close supervisions is necessary when this appliance is used by on, or near children, invalids or disabled persons.
- Use this appliance only for its intended use as described in this manual. Do not use any attachments not provided with the product.
- **NEVER** operate the appliance if it has a damaged cord or plug, if it is not working properly, been dropped or damaged, or dropped into water, contact our After Sales Support.
- **NEVER** drop or insert any object into any opening.
- **NEVER** operate the appliance with the air opening blocked. Keep the air openings free of lint and hair.
- **DO NOT** operate where aerosol (spray) products are being used or where oxygen is being administered.
- **DO NOT** operate under a blanket or pillow. Excessive heating can occur and cause fire, electric shock or injury to persons.
- **DO NOT** carry this appliance by supply cord or use the cord as a handle.
- **DO NOT** use outdoors.
- **DO NOT** stand on or in appliance. Use only when seated.
- **DO NOT** handle the plug with wet hands
- This product should never be used by any individual suffering from any physical alignment that would limit the capacity to operate the controls, or who have sensory deficiencies.

## Important Safety Instructions (Cont.)



- **ALWAYS** unplug the unit when filling or moving the unit.
- Ensure that your hands are dry when operating the controls or removing the plug.
- Use the appliance on a sturdy or level surface.
- This appliance has a heated surface. Persons insensitive to heat must be careful when using the appliance.
- This product is for indoor household use only.

### Health matters

- This product is not intended for medical use.
- Consult your doctor prior to using this product if:
  1. You are pregnant
  2. You have a pacemaker
  3. You have concerns regarding your health
- **NOT** recommended for use by diabetics.
- **NEVER** use the device on any open wounds, discoloured areas, or any area of the body that is swollen, burned or inflamed.
- **NEVER** use the product after surgery or if you are not yet fully recovered.
- If you experience pain in a muscle or joint for an extended period of time, discontinue use and consult your physician.

- **Do not** use the product in a high-temperature environment or where it may become wet or damp. If the unit becomes wet, do not use.
- **Do not** use this product if you are currently pregnant or menstruating.
- **Do not** place the unit in a place where it may be in contact with flammable items, explosive items, or chemicals.
- **Do not** use this product if, you have any medical complication such as osteoporosis, a skin injury or disorder, low blood pressure, cancer or a pacemaker or other medical implant. Contact your physician before using this product.
- Stop using the product immediately if you feel uncomfortable or any abnormalities occur.

# Product Overview

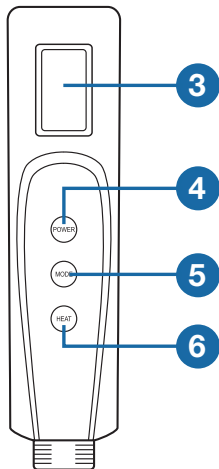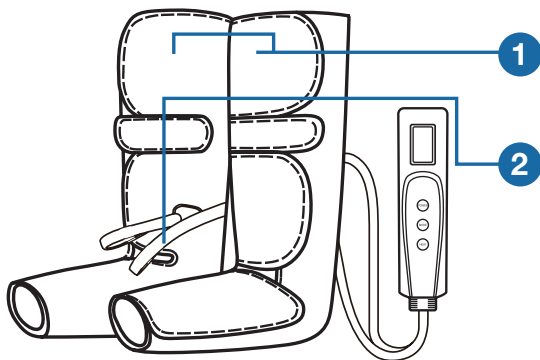

- 1 Leg massager
- 2 Air pressure hose
- 3 LCD screen
- 4 Power button: press and hold to turn the device on. Press and hold again to turn the device off. Press repeatedly to change between intensity modes (low/middle/high intensity mode).
- 5 Mode button: press to choose between the 3 built-in modes.
- 6 Heat button: press to turn heating function on/off (note: heating will only work on bottom of feet).

# Operation

- Slip your leg onto the leg massagers, keeping your feet flat on the bottom. Each leg massager is marked L and R.
- Secure your leg into each massager by tightening the strap.
- Ensure that the massager is securely fastened on both legs.
- Connect the air pressure hose to the remote control.
- Press and hold the Power button on the remote control to turn the device on. The LCD screen will light up.
- Select the desired intensity and mode with the Power and Mode buttons on the remote control.
- Mode 1: Calf + Foot airbag massage
- Mode 2: Calf massage
- Mode 3: Foot airbag massage
- **Note:** the device will automatically shut off after 20 minutes.
- **Note:** it is recommended to use the device no longer than an hour.

# Troubleshooting

Do not attempt to open and repair the product yourself.

Problem	Solution
Device is not working	Ensure the device is connected to a power source. Press and hold the Power button on the remote control. The LCD screen will light up.
Device stopped working during massage	The device will automatically shut off after 20 minutes. Restart the device and resume use.



# Specifications

Power Adaptor input	AC 100-240V, 50/60Hz, 0.5A Max.
Power Adaptor output	12V, 1.5A
Power consumption	18W
Weight	1.2 kg (approx.)

For Technical and Warranty Support

1300 886 649

[tempo.org/support](http://tempo.org/support)

Distributed by Tempo (Aust) Pty Ltd,  
PO Box 132, Frenchs Forest NSW 1640

## Compliance

This appliance has been fully tested and meets all requirements as set out by standards AS/NZS 60335.2.32 for electrical safety.



The RCM Mark (Regulatory Compliance Mark) indicates that the product complies with the relevant guidelines of the ACMA as well as corresponding government requirements for the safety of electrical devices.

# Notes

# Notes



## **Warranty returns**

Should you for any reason need to return this product for a warranty claim, make sure to include all accessories with the product.

## **Product does not work?**

If you encounter problems with this product, or if it fails to perform to your expectations, make sure to contact our After Sales Support Centre on (AU) 1300 886 649 or (NZ) 0800 836 761 for advice.

Distributed by Tempo (Aust) Pty Ltd ABN 70 106 100 252

PO BOX 132 Frenchs Forest, Australia NSW 1640

Customer Helpline:

(AU) 1300 886 649 (NZ) 0800 836 761

Email: [tempo.org/support](mailto:tempo.org/support)